




MGF 2 1ÖS/1ÖS-MED

Part-No.: 1441015

Specification

- Footswitch in light grey RAL7035
- 2 Actuators (B7F) in signal blue RAL5005
- Cable in light grey RAL7035
- With bracket

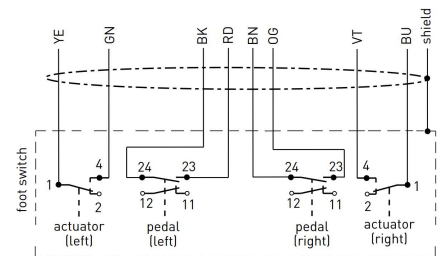

Technical data

Standards	IEC 60601-1 (Basic Safety) IEC 60529 (Degrees of protection) Note: see CSA certificate for corresponding revision level.
Enclosure + cover	GD-Al alloy, coloured in light grey RAL7035
Pedal	GD-Al alloy, coloured in light grey RAL7035
Actuator B7F	RAL5005 signal blue
Bracket	electro polished
Degree of protection	IPX8 (1m / 35 min.) as per DIN EN 60529
Contact element and material	1NC/1NO per pedal (wired as normally open) / gold-plated
Switch insert	Slow action, forced contact
Actuator B7F Contact element and material	Micro Switch BURGESS V4NCBT7 / Silver
Connection	3 m shielded cable (8xAWG24 C UL, RAL7035), without connector
Mechanical lifetime	The expected operating life of the foot switch is limited to 10 years or 1,000,000 switching cycles of the actuating units, whichever comes first.
Dimensions (W x L xH)	195 x 221 x 61 mm (without bracket) 224 x 221 x 171 mm (with bracket)
Weight	0,860 kg



MGF 2 1ÖS/1ÖS-MED

Part-No.: 1441015

Electrical ratings actuator

Switching voltage (U_B)	25 VAC / 60VDC
Switching current (I_B)	100 mA
Switching capacity (P_B)	6 W

Environmental conditions

	<u>Storage / Transport</u>	<u>Operating (Operation)</u>
Temperatur range	-40 °C up to +70 °C	-10 °C up to +60 °C
Relative humidity	10 % up to 100 %	10 % up to 100 %
Air pressure	500 hPa – 1120 hPa	800 hPa – 1060 hPa

Mounting/Wiring

Wire the foot switches according to the specified wire colours/terminal labelling. Contact symbols are shown for a not actuated switch.

Maintenance

Depending on the environmental/operating conditions we recommend routine maintenance as follows:

- Check enclosure and connecting cable for damage and destructive dirt
- Check functional elements and actuating elements for free operation
- For only manual cleaning, use just a cloth impregnated with water and a mild detergent.
- Do not use cleaning agents which may damage the plastic surfaces such as detergents, abrasive cleansers or solvent-based cleaners (such as benzine, stain remover).

Disposal

- Observe national, local and legal regulations concerning disposal.
- Recycle each material separately.

General notes

The foot switch as a component of a medical device can be evaluated only together in the overall system of the customer. Therefore the conformity evaluation including the classification of the overall system according to the Medical Device Directive 93/42/EEC respectively EU Medical Device Regulation 2017/745 must be carried out at the customer. In addition, validation can also only be carried out by the customer in the overall system.

Furthermore, the intended use and the basic functions are verifiable only in interaction with the overall system. As a result, the customer must involve the footswitch in his risk management, his usability analysis, the verification of the electromagnetic compatibility and, if applicable, the verification of the biocompatibility.



MGF 2 1ÖS/1ÖS-MED

Part-No.: 1441015

Safety instructions

Despite careful production and inspection as well as redundant design of important components, an entirely safe function cannot be guaranteed at all times (in case of malfunction, for example, a switch contact can turn off incorrectly, or a spring may break).

A functional test must be performed before every operation.

Risk of injury due to incorrect handling!

- The plug connections must not be subjected to mechanical loads!
- Switch off terminal before cleaning or maintenance work!
- Test functions when device is switched back on!
- Moving actuating components (e.g. carrying handle) can be a crushing hazard!
- Tripping hazard!

If components with reed contacts or hall sensors are used, strong magnetic fields can result in unwanted influences to the system. Additionally, magnetically conducting materials are to be kept away from the direct foot switch area!

The electrical connection may only be carried out by authorized personnel.

Components, especially pedals or other actuating elements, may not be disassembled or dismembered under any circumstances! If this happens the products need to be sent back immediately.

Errors and omissions excepted.